

APA222Mu01 50µg

Active Interferon Beta (IFNb)

Organism Species: Mus musculus (Mouse)

Instruction manual

FOR RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

12th Edition (Revised in Aug, 2016)



[PROPERTIES]

Source: Prokaryotic expression

Host: E.coli

Residues: Ile22~Asn182

Tags: N-terminal His Tag

Purity: > 95%

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 1mM DTT, 5% Trehalose and

Proclin300.

Applications: Cell culture; Activity Assays.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.7

Predicted Molecular Mass: 49.7kDa

Accurate Molecular Mass: 44&48kDa as determined by SDS-PAGE reducing conditions.

Phenomenon explanation:

The possible reasons that the actual band size differs from the predicted are as follows:

- 1. Splice variants: Alternative splicing may create different sized proteins from the same gene.
- 2. Relative charge: The composition of amino acids may affects the charge of the protein.
- 3. Post-translational modification: Phosphorylation, glycosylation, methylation etc.
- 4. Post-translation cleavage: Many proteins are synthesized as pro-proteins, and then cleaved to give the active form.
- 5. Polymerization of the target protein: Dimerization, multimerization etc.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.


Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.



Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

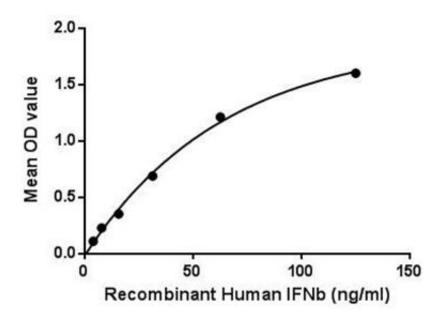

[IDENTIFICATION]

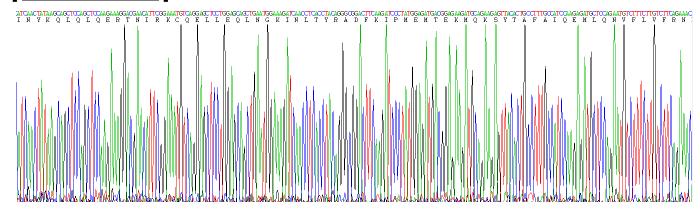


Figure . Gene Sequencing (extract)



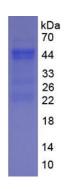


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.