


MAC979Hu22
Monoclonal Antibody
To Carcinoembryonic Antigen Related Cell Adhesion Molecule 7 (CEACAM7)
Organism Species: Homo sapiens (Human)
Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

9th Edition (Revised in Jul, 2013)

[PRODUCT INFORMATION]

Immunogen: CEACAM7, Human
Clonality: Monoclonal
Clone number: G1
Host: Mouse
Immunoglobulin Type: IgG
Purification: Affinity Chromatography.
Applications: WB, ICC, IHC-P, IHC-F, ELISA
Concentration: 500µg/mL
UOM: 200µg

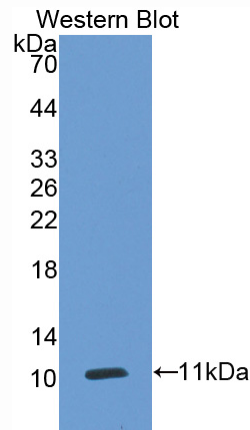

Sample: Recombinant CEACAM7, Human

[IMMUNOGEN INFORMATION]

Immunogen: Recombinant CEACAM7 (Lys147~Thr231) expressed in *E.coli*.
Accession No.: RPC979Hu01
Sequence: The target protein is fused with N-terminal His-Tag and its sequence is listed below.

MGHHHHHSGSEF-KPSI TSNNFNPVEN KDIVLTCQP ETQNTTYLWW VNNQSLLVSP
RLLLSTDNRT LVLLSATKND IGPYECEIQN PVGASRSDPV T

[ANTIBODY SPECIFICITY]

The antibody is a mouse monoclonal antibody raised against CEACAM7. It has been selected for its ability to recognize CEACAM7 in immunohistochemical staining and western blotting.

[APPLICATIONS]

Western blotting: 1:100-400

Immunocytochemistry in formalin fixed cells: 1:100-500

Immunohistochemistry in formalin fixed frozen section: 1:100-500

Immunohistochemistry in paraffin section: 1:50-200

Enzyme-linked Immunosorbent Assay: 1:100-200

Optimal working dilutions must be determined by end user.

[CONTENTS]

Form & Buffer: Supplied as solution form in PBS, pH7.4, containing 0.02% NaN₃, 50% glycerol.

[STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.