

**rP90815Ra01 10 $\mu$ g**  
**Nitric Oxide Synthase 1, Neuronal (NOS1)**  
**Organism: Rattus norvegicus (Rat)**  
***Instruction manual***

FOR IN VITRO USE AND RESEARCH USE ONLY  
NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

7th Edition (Revised in May, 2013)

## [ **PROPERTIES** ]

**Residues:** Gln9~Ser136 (Accession # P29476),  
with N-terminal His-Tag.

**Host:** *E. coli*

**Subcellular Location:** Cell membrane, dendritic  
spine.sarcolemma; Peripheral membrane protein.  
Cell projection.

**Purity:** >95%

**Endotoxin Level:** <1.0EU per 1 $\mu$ g  
(determined by the LAL method).


**Formulation:** Supplied as lyophilized form in PBS,  
pH7.4, containing 5% sucrose, 0.01% sarcosyl.

**Predicted isoelectric point:** 7.1

**Predicted Molecular Mass:** 15.2kDa

**Applications:** SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

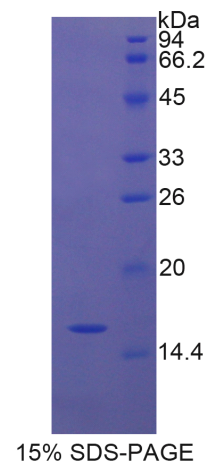

## [ **USAGE** ]

Reconstitute in sterile PBS, pH7.2-pH7.4.

*Unique product Superb quality Client favorite Nicest service*  ISO9001:2008;  ISO13485:2003; 

## [ STORAGE AND STABILITY ]

### **Storage: Avoid repeated freeze/thaw cycles.**

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

**Stability Test:** The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

## [ SEQUENCES ]

The target protein is fused with N-terminal His-Tag, its sequence is listed below.  
MGHHHHHSGSEF-QQ IQPNVISVRL FKRKVGGLGF LVKERVSKPP VIISDLIRGG  
AAEQSGLIQA GDIILAVNDR PLVDLSYDSA LEVLRGIASE THVVLILRGP EGFTTHLETT  
FTGDGTPKTI RVTQPLGPPT KAVDLS

## [ REFERENCES ]

1. Hugyec M., *et al.* (2011) Brain Res. 1404:31-38.
2. Hyndman K.A., *et al.* (2011) Am. J. Physiol. Renal Physiol. 301:F118-24.
3. Talapka P., *et al.* (2011) World J. Gastroenterol. 17:1026-1029.
4. Aoki T., *et al.* (2011) Lab. Invest. 91:619-626.