

PAB765Mu01 Polyclonal Antibody To Killer Cell Lectin Like Receptor Subfamily K, Member 1 (KLRK1) Organism Species: Mus musculus (Mouse) Instruction manual FOR IN VITRO USE AND RESEARCH USE ONLY

FOR IN VITRO USE AND RESEARCH USE ONLY NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

### [PRODUCT INFORMATION]

Immunogen: KLRK1, Mouse Clonality: Polyclonal Host: Rabbit Immunoglobulin Type: IgG Purification: Affinity Chromatography. Applications: WB, ICC, IHC-P, IHC-F, ELISA Concentration: 200µg/mL UOM: 100µg

#### [IMMUNOGEN INFORMATION]

**Immunogen:** Recombinant KLRK1 (Phe90~Val232) with two N-terminal Tags, His-tag and T7-tag expressed in *E.coli*.

Accession No.: RPB765Mu01

### [ANTIBODY SPECIFITY]

The antibody is a rabbit polyclonal antibody raised against KLRK1. It has been selected for its ability to recognize KLRK1 in immunohistochemical staining and western blotting.

## [APPLICATIONS]

Western blotting: 1:50-400 Immunocytochemistry in formalin fixed cells: 1:50-500 Immunohistochemistry in formalin fixed frozen section: 1:50-500 Immunohistochemistry in paraffin section: 1:10-100 Enzyme-linked Immunosorbent Assay: 1:100-200 Optimal working dilutions must be determined by end user.

# Coud-Clone Corp.

## [CONTENTS]

**Form & Buffer:** Supplied as solution form in PBS, pH7.4, containing 0.02% NaN<sub>3</sub>, 50% glycerol.

## [QUALITY CONTROL]

**Content:** The quality control contains recombinant KLRK1 (Phe90~Val232) disposed in loading buffer.

**Usage:** 10uL per well when 3,3'-Diaminobenzidine(DAB) as the substrate.


5uL per well when used in enhanced chemilumescent (ECL). **Note:** The quality control is specifically manufactured as the positive control. Not used for other purposes.

**Loading Buffer:** 100mM Tris(pH8.8), 2% SDS, 200mM NaCl, 50% glycerol, BPB 0.01%, NaN<sub>3</sub> 0.02%.

## [STORAGE]

Store at 4°C for frequent use. Stored at -20°C to -80°C in a manual defrost freezer for one year without detectable loss of activity. Avoid repeated freeze-thaw cycles.

## [<u>IMAGES</u>]

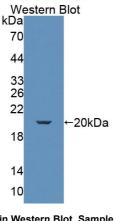

Used in DAB staining on fromalin fixed paraffin- embedded spleen tissue

Used in Western Blot, Sample: Recombinant KLRK1, Mouse

1304 Langham Creek Dr, Suite 226, Houston, TX 77084, USA | 001-888-960-7402 | www.cloud-clone.us | mail@cloud-clone.as Export Processing Zone, Wuhan, Hubei 430056, PRC | 0086-800-880-0687 | www.cloud-clone.com | mail@cloud-clone.com