

RPF466Ra01 100µg

Recombinant Centaurin Alpha 2 (CENTa2)

Organism Species: *Rattus norvegicus* (Rat)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Ser49~Ser360

Tags: N-terminal His Tag

Subcellular Location: Membrane, Cytoplasm

Purity: > 80%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose .

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 9.6

Predicted Molecular Mass: 40.4kDa

Accurate Molecular Mass: 42kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

| | | | | | |
|------------|------------|------------|------------|------------|----|
| VHRNFPDISK | VKSVRLDFWD | DSMVEFMTHN | GNLSVKAKFE | ARVPTFYYPV | SG |
| QASDCLVLKE | QWIRAKYERQ | EFMAEKAVSP | PGDREGFLWK | RGRDNSQFLR | |
| RRFVLLSREG | LLKYTTKEEG | KTPKAIINIK | DLNATFQTEK | IGHPHGLQIT | |
| YRKEGQTRNL | FVYHDSGKEI | VDWFNALRAA | RLQYLKLAFF | DLPESELVPL | |
| ITRNYLKQGF | MEKTGPKHRE | PFKKRWFALD | PQERRLLYYK | NPLDAFEQQQ | |
| VFLGSNEQGY | EVWEGLPQGI | RGNRWKVGLT | VITPERKFVF | TCPTEKEQRE | |
| WLESLRGVLS | | | | | |

[IDENTIFICATION]

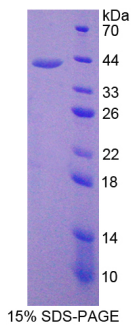

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.