

RPA800Ra01 100µg

Recombinant Pregnancy Associated Plasma Protein A (PAPPA)

Organism Species: *Rattus norvegicus* (Rat)

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Ala1111~Gly1440

Tags: N-terminal His Tag

Subcellular Location: Secreted

Purity: > 90%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% SKL, 5% Trehalose.

Original Concentration: 200µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 6.8

Predicted Molecular Mass: 40.3kDa

Accurate Molecular Mass: 40kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in 10mM PBS (pH7.4) to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.


Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

ACEPVDCGIP DHHHVYAASF SCPEGNNSFL TCMEDGLWSF PEALCELMCL
APPPVFNADL QTARCRENKH KVGSLCKYKC KPGYHVPSS RKSCKRAFKT
QCTQDGSWQE GTCVPVTCDF PPPKFHGLYQ CTNGFQFNSE CRIKCEDSDA
SQGRGSNTIH CRKDGTWSGS FHVCREMQGQ CSAPNQLNSH LKLQCPDGYA
IGSECATSCL DHNSESILP VNLTVRDIPH WMNPTRVQRI VCTAGLQWYP
HPARIHCVKG CEPFMGDNYC DAINNRAFCN YDGGDCCTST VKTKKVTPFP
MSCDLQDDCA CRDPEAQEHS RKDLRGYSHG

[IDENTIFICATION]

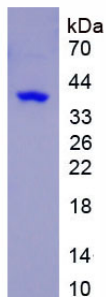

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was used in clinical diagnostic or any other procedures.