

RPL074Hu01 100µg

Recombinant Ubiquitin Protein Ligase E3 Component N-Recognin 4 (UBR4)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Val124~Arg396

Tags: Two N-terminal Tags, His-tag and T7-tag

Accession: Q5T4S7

Host: *E. coli*

Subcellular Location: Membrane.

Purity: >95%

Endotoxin Level: <1.0EU per 1µg

(determined by the LAL method).


Formulation: Supplied as lyophilized form in PBS, pH7.4, containing 5% trehalose, 0.01% sarcosyl.

Predicted isoelectric point: 5.3

Predicted Molecular Mass: 33.5kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

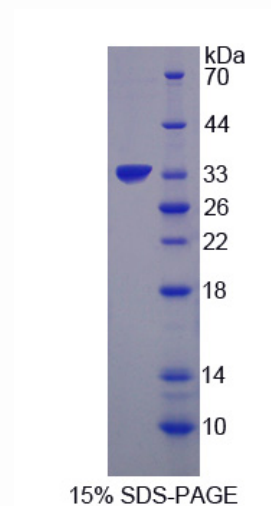

[USAGE]

Reconstitute in sterile PBS, pH7.2-pH7.4.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

VSQKHLI LLIKGLCTGC SRLDRTEIIT FTAMMKSACL PQTVKTLSDV EDQKELASPV
SPELRQKEVQ MNFLNQLTSV FNPRTVASQP ISTQTLVEGE NDEQSSTDQA SAIKTKNVFI
AQNVASLQEL GGSEKLLRVC LNLPHYFLRYI NRFQDAVLAN SFFIMPATVA DATAVRNGFH
SLVIDVTMAL DTLSLPVLPEP LNPSRLQDVT VLSLSCLYAG VSVATCMAIL HVGSAQQVRT
GSTSSKEDDY ESDAATIVQK CLEIYDMIGQ AISSSR