

RPC099Hu01 100µg

Recombinant Integrin Beta 6 (ITGb6)

Organism Species: Homo sapiens (Human)

Instruction manual

FOR IN VITRO USE AND RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

10th Edition (Revised in Jan, 2014)

[PROPERTIES]

Residues: Tyr131~Leu371

Tags: Two N-terminal Tags, His-tag and T7-tag

Accession: P18564

Host: *E. coli*

Subcellular Location: Membrane; Single-pass type I membrane protein.

Purity: >95%

Endotoxin Level: <1.0EU per 1µg
(determined by the LAL method).

Formulation: Supplied as lyophilized form in 20mM Tris, 150mM NaCl, pH8.0, containing 0.01% sarcosyl and 5% trehalose.

Predicted isoelectric point: 4.8


Predicted Molecular Mass: 30.5kDa

Applications: SDS-PAGE; WB; ELISA; IP.

(May be suitable for use in other assays to be determined by the end user.)

[USAGE]

Reconstitute in sterile ddH₂O.

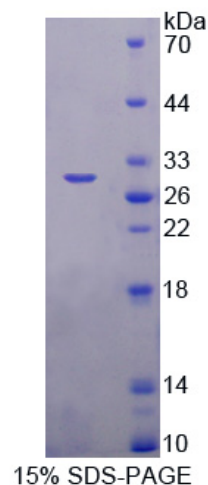

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate of the target protein. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. (Referring from China Biological Products Standard, which was calculated by the Arrhenius equation.) The loss of this protein is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCES]

The sequence of the target protein is listed below.

YPVDLYLMD LSASMDDDLN TIKELGSRLS KEMSKLTSNF RLGFGSFVEK PVSPFVKTTTP
EEIANPCSSI PYFCLPTFGF KHILPLTND AERFNEIVKNQ KISANIDTPE GGFDAIMQAA
VCKEKIGWRN DSLHLLVFVS DADSHFGMDS KLAGIVIPND GLCHLDSKNE YSMSTVLEYP
TIGQLIDKLV QNNVLLIFAV TQEQVHLYEN YAKLIPGATV GLLQKDSGNI LQLIISAYEE L

[REFERENCES]

1. Sheppard D., *et al.* (1990) J. Biol. Chem. 265:11502-11507.
2. Jiang W.-M., *et al.* (1992) Int. Immunol. 4:1031-1040.
3. Agrez M.V., *et al.* (1997) Virology 239:71-77.
4. van Der Flier A., *et al.* (2002) J. Cell Biol. 156:361-376.