

RPT479Hu01 50µg

Recombinant Tripartite Motif Containing Protein 56 (TRIM56)

Organism Species: *Homo sapiens (Human)*

Instruction manual

FOR RESEARCH USE ONLY

NOT FOR USE IN CLINICAL DIAGNOSTIC PROCEDURES

13th Edition (Revised in Aug, 2023)

[PROPERTIES]

Source: Prokaryotic expression

Host: *E.coli*

Residues: Met1~Ser755

Tags: N-terminal His Tag

Subcellular Location: Nucleus, Cytoplasm

Purity: > 80%

Traits: Freeze-dried powder

Buffer formulation: PBS, pH7.4, containing 0.01% Sarcosyl, 5% Trehalose.

Original Concentration: 100µg/mL

Applications: Positive Control; Immunogen; SDS-PAGE; WB.

(May be suitable for use in other assays to be determined by the end user.)

Predicted isoelectric point: 7.8

Predicted Molecular Mass: 85.2kDa

Accurate Molecular Mass: 85kDa as determined by SDS-PAGE reducing conditions.

[USAGE]

Reconstitute in ddH₂O to a concentration of 0.1-1.0 mg/mL. Do not vortex.

[STORAGE AND STABILITY]

Storage: Avoid repeated freeze/thaw cycles.

Store at 2-8°C for one month.

Aliquot and store at -80°C for 12 months.

Stability Test: The thermal stability is described by the loss rate. The loss rate was determined by accelerated thermal degradation test, that is, incubate the protein at 37°C for 48h, and no obvious degradation and precipitation were observed. The loss rate is less than 5% within the expiration date under appropriate storage condition.

[SEQUENCE]

MVSHGSSPSLLEALSSDFLACKICLEQLRAPKTLPC LH TYCQDCLAQLADGGRVRCPECRET
 VPVPPEGVASFKTNFFVNGLLDLVKARACGDLRAGKPACALCPLVGGTSTGGPATARCLDCA
 DDLCQACADGHRCTRQTHTHRVDLVGYRAGWYDEEARERQAAQCPQHPGEALRFLCQ
 PCSQLLCRECRDPHLDHPCLPLAEAVRARRPGLEGLLAGVDNNLVELEAARRVEKEALARL
 REQAARVGTQVEEAAEGVLRALLAQKQEV LGQLRAHVEAAEEAARERLAELEGREQVARA
 AAFAARRVLSL GREAEILSLEGAIAQRLRQLQGCPWAPGPAPCLLPQLELHPGLLDKNCHLL
 RLSFEEQQPQKDGGKDGAGTQGGEESSRREDEPKTERQGGVQPQAGDGAQTPKEEKA
 QTTREEGAQTLEEDRAQTPHEDGGPQPHRGGRPNKKKKFKGRLKSSREPSALGPNLDGS
 GLLPRPIFYCSFPTRM PGDKRSPRITGLCPFGPREILVADEQNRALKRFSLNGDYKGTVPVPEG
 CSPCSVAALQSAVAFSASARLYLINPNGEVQWRRALSLSQASHAVAALPSGDRVAVSVAGHV
 EVYNMEGSLATRFIPGGKASRGLRALVLTTS PQGHFVGSDWQQNSVVICDGLGQVVGEYK
 GPGLHGCQPGSVSVDKKGYIFLTLEVNKVVILDPKGSLLGDFLTAYHGLEKPRVTTMVDGRY
 LVVLSLNGTIHIFRVRSPDS

[IDENTIFICATION]

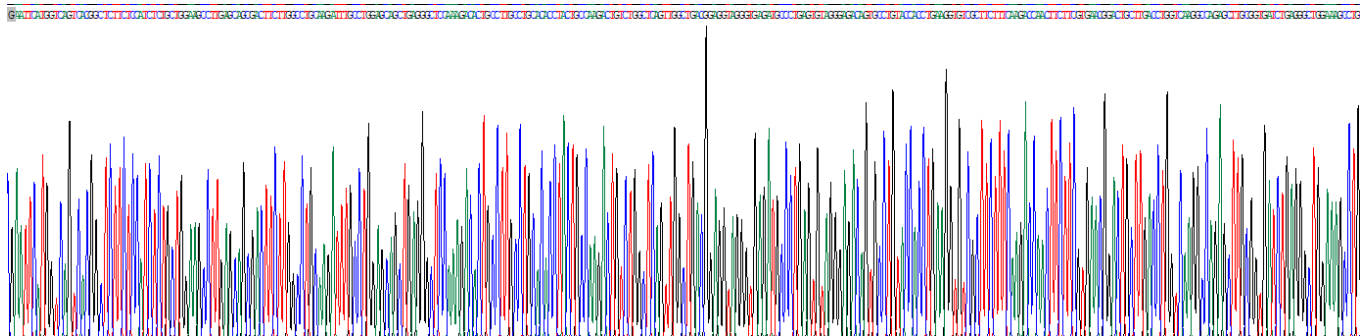


Figure . Gene Sequencing (extract)

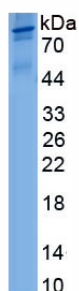


Figure. SDS-PAGE

[IMPORTANT NOTE]

The kit is designed for research use only, we will not be responsible for any issue if the kit was

used in clinical diagnostic or any other procedures.